
niceboy ION[®]

ION SONIC PRO UV

User Manual / Sonic Brush

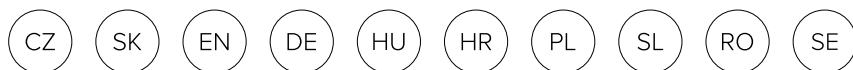OBSAH BALENÍ

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

PACKAGE CONTENTS

- Toothbrush handle
- 3x brush head (soft, medium, hard)
- Magnetic docking station
- USB-C charging cable
- Rechargeable travel case with UV function
- Instructions for use

SCHEMATIC OF THE SONIC UV PRO PRODUCT

1. Bristle
2. The header
3. Body
4. On/off/mode switch button
5. Indicator of cleaning modes
6. Power on/charge indicator
7. Charging connector

CLEANING MODES

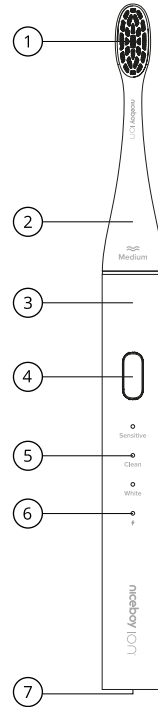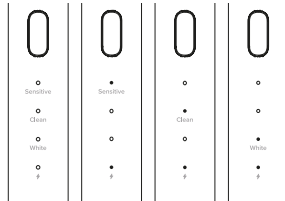INSTRUCTIONS FOR USE


1. When using for the first time, place the brush in the magnetic docking station.

Charging process:


- 0%-10%, charging indicator shows steady red color
- 10-100%, the charging indicator shows a steady orange color
- Charging is complete, the charging indicator shows steady green

Note: If the product has not been used for a long time, the charging indicator may not appear when charging the electric toothbrush. Please wait patiently for 3-5 minutes to resume normal charging.

2. Insert the new brush head as shown; if you need to take it out, just pull it out as shown in the picture.
3. Wet the bristles and squeeze out an appropriate amount of toothpaste. Touch the bristles coated with toothpaste to the surface of the tooth in any part of the mouth.
4. Press the switch and move the jaw slightly up and down to clean your teeth laterally from the outside as well as from the inside. Slowly move the brush bristles at a 45° angle.


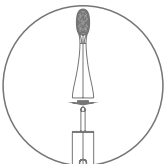
5. The oral cavity is divided into 4 areas and the recommended cleaning time for each area is 30 seconds. After you finish cleaning that area for 30 seconds, the electric toothbrush will pause briefly to remind you to clean the next area .

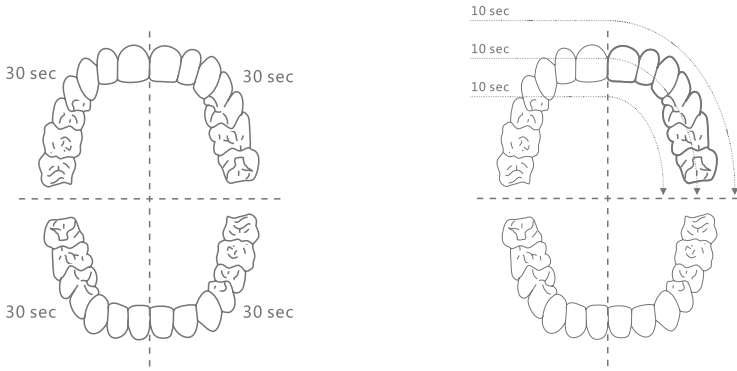6. When cleaning, do not exert unnecessary pressure on the teeth and do not make circular movements of the brush. Cleaning with a sonic toothbrush should be done with slow movements over the surface of the teeth.

7. After the scientifically determined brushing time of 2 minutes, the electric toothbrush will automatically stop. After cleaning the teeth, rinse the brush head with running water, shake off excess water from the brush head and handle. Keep toothbrush crystals dry and store in a dry place.

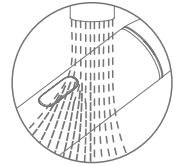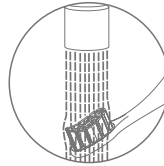
Note:

- If you are not going to use the brush for a long time, disconnect the charging connector and thoroughly clean the head and handle of the brush, store in a cool and dry place;
- The head replacement interval is recommended after 60 days.
- Auto memory mode: use any mode for more than 5 seconds, the system will record the mode. By default, memory mode is activated the next time you turn it on.

TRAVEL CASE

1. Charging indicator
2. Charging USB-C connector
3. UVC on indicator

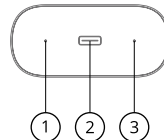
TRAVEL CASE INFORMATION

- To charge and activate UVC, it is necessary to connect the case to the network using the USB-C connector.
- It takes 30 seconds to sterilize the brush head.


INSTRUCTIONS FOR USING THE TRAVEL CASE


1. Open the case by sliding the top side in the direction of the arrow in the following image.

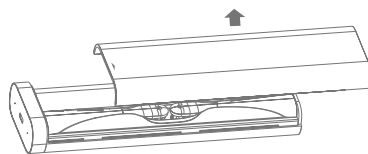2. Remove the top of the case.

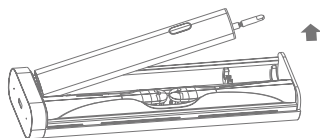3. Insert/remove the brush body from/into the case as shown in the following image.

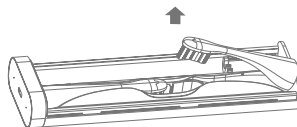4. Insert/remove the brush heads from/into the case as shown in the following figure.

ADDITIONAL INFORMATION ABOUT THE TRAVEL CASE

- To achieve 99.9% sterilization of the brush, it is first necessary to let the head dry and then put it in the UVC case.

SAFETY INSTRUCTIONS

- Please read this instruction manual carefully before using the appliance and keep it for future reference.
- Keep the charging base plugged in away from water. Do not place or store it near or near water contained in a bathtub, sink, etc. Do not immerse the charging base in water or any other liquid. After cleaning, make sure the charging base is completely dry before connecting it to the mains.
- Before connecting the appliance, check that the voltage indicated on the appliance corresponds to the voltage in the local electricity network.
- Never use the charging base if it is damaged in any way. If the charging base cable or connector is damaged, always have it replaced with a new one to avoid danger.
- If the appliance is damaged in any way (brush head, toothbrush handle), stop using it. This appliance does not contain any user-serviceable parts. If the appliance is damaged, contact the consumer care center in your country.
- Do not use the charging base outdoors or near heated surfaces.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.

5. PRELIMINARY MEASURES

- Do not use the toothbrush in the bath or shower.
- If you have had gum or oral surgery in the previous 2 months, consult your dentist before using the toothbrush.
- If excessive bleeding occurs after using this toothbrush, or if bleeding continues after 1 week of use, consult your dentist.
- The toothbrush meets safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, check with your doctor or the manufacturer of the implanted device before using the brush.
- This appliance was designed only for cleaning the teeth, gums and tongue. Do not use it for any other purpose. If you experience any discomfort or pain, stop using the appliance and contact your doctor.
- If you have health problems, consult your doctor before using the device. - The toothbrush is a personal care device and is not intended for use by multiple patients in a dental office or facility.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear. Do not clean the brush head in the dishwasher.
- Do not use brush heads other than those recommended by the manufacturer.
- Supervise children while brushing until they develop proper brushing habits.
- If the toothpaste contains peroxide, baking soda, or bicarbonate (common in whitening toothpastes), clean the brush head thoroughly with soap and water after each use. This will prevent possible cracking of the plastic.

SONIC PRO UV PRODUCT SPECIFICATIONS

RPM:	85,000
Voltage / power:	3.7V ~ 1.85W
Charging:	Magnetic charger with USB-C input
Input voltage:	5V ~ 300 mA
Operating time:	55 days (calculated according to the operating status of 2 minutes/time and 2 times/day)
Battery capacity:	800 mAh
Charging time:	3 hours
Waterproof:	IPX7
Recommended charger:	USB adapter with a voltage of 5V/1A or 5V/2A

SPECIFICATIONS OF THE PRODUCT SONIC PRO UV TRAVEL CASE

Voltage / power:	3.7V ~ 0.4W
Charging:	USB-C cable
Input voltage:	5V ~ 500mA
Recommended charger:	USB adapter with a voltage of 5V/1A or 5V/2A

MORE INFORMATION

The company Niceboy s.r.o. hereby declares that the Sonic Pro UV radio equipment type is in compliance with directives 2014/30/EU, 2014/35/EU and 2011/65/EU. The full content of the EU Declaration of Conformity is available on the following websites:
<https://niceboy.eu/en/declaration/ion-sonic-pro-uv>

INFORMATION FOR USERS FOR THE DISPOSAL OF ELECTRICAL AND ELECTRONIC EQUIPMENT (HOMES)

 This symbol placed on the product or in the original product documentation means that used electrical or electronic products must not be disposed of together with municipal waste. If you want to dispose of these products correctly, take them to a designated collection point, where they will be accepted free of charge. By disposing of the product in this way, you help protect precious natural resources and help prevent possible negative impacts on the environment and human health that could result from improper waste disposal. Please contact your local authority or the nearest collection point for more detailed information. According to national regulations, fines can also be issued to anyone who disposes of this type of waste incorrectly. Information for users regarding the disposal of electrical and electronic equipment (Business and Corporate Use). For proper disposal of electrical and electronic equipment for commercial and business use, contact the manufacturer or importer of the product. They will provide you with information on all disposal methods and, according to the date indicated on the electrical or electronic equipment on the market, tell you who is responsible for financing the disposal of that electrical or electronic equipment. Information regarding disposal procedures in other non-EU countries. The above symbol only applies to European Union countries. For the correct disposal of electrical and electronic equipment, ask your local authorities or the equipment dealer for the appropriate information.

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